
ingo-man® smart

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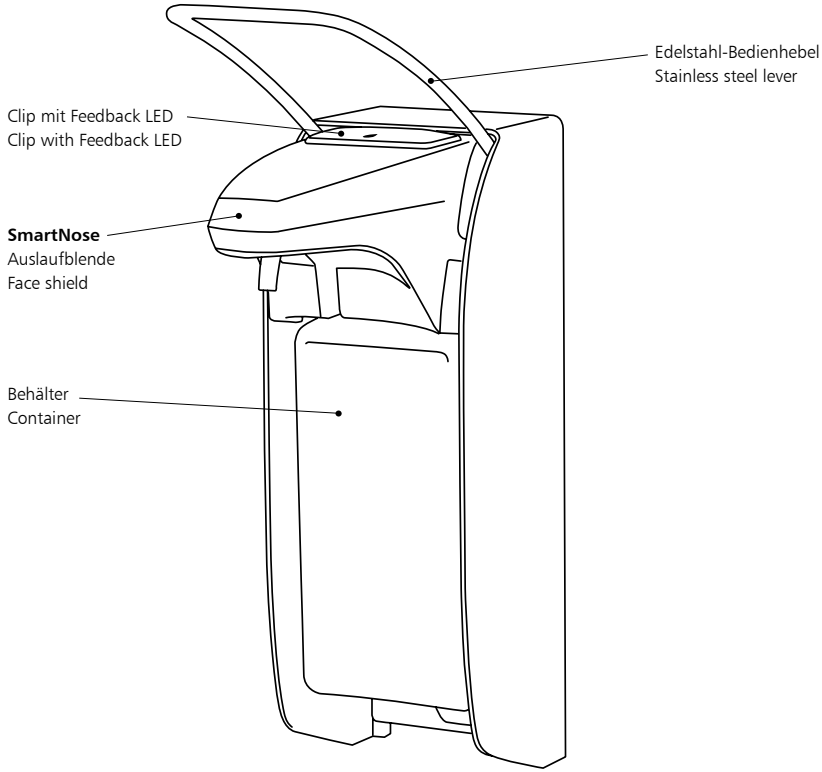E General safety

Repairs to the appliance may exclusively be carried out by OPHARDT Hygiene employees or persons authorized by OPHARDT Hygiene. Important dangers may result from by unexpert repairs. In case of a necessary repair please contact the manufacturer.

Locate the appliance in a dry and splash-proof place.

This appliance is not a toy. Keep and use it out of the reach of children. The appliance should be installed where children cannot reach it.

Operate the device only with an original battery, which is also available as a spare part. (Lithium battery 9 Ah 3.6 V, Art. No. 4402691 or Lithium battery 3.5 Ah 3.6 V, Art. No. 3401863. See page 27 for details.)

This device contains a non-rechargeable lithium battery.

Do not operate the appliance with rechargeable batteries (accumulators).

Never use force.

Dispose of the empty batteries in line with local policy and protocol.

The device may only be opened to change the battery.

Do not expose the device to thermal radiation or heat.

E Safety during operation of the appliance

Protect the appliance against humidity, the device must not be sprayed for disinfection/ cleaning.

Never submerge the appliance in water.

There are no dismountable parts with the exception of SmartNose, magnet holder, battery.

The appliance is not autoclavable and cannot be cleaned in the dishwasher.

The device must not be cleaned with water jets or by means of a high-pressure cleaner.

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Intended use

The intelligent face shield *SmartNose* records hand hygiene behaviour fully automatically and continuously.

The device is exclusively designed for battery operation.

The device may only be used indoors.

Ambient conditions:

- Never use the appliance in potentially explosive areas.
- Only use the appliance in a dustless environment.
- Pay attention that the apparatus is only used at a relative air moisture of <90% (not condensing).
- Only use this apparatus at room temperature

Operation of the SmartNose according to the regulations

In order to always guarantee optimum functioning of the dispenser the following must be observed:


- When replacing a battery please absolutely pay attention to the correct polarity! (Reverse polarity protected plug connection)
- The SmartNose may not be located close to radiation-emitting or statically charged objects.

Attention!


OPHARDT does not assume any liability for damages in the case of:

- defects to the appliance caused by mechanical influences or overvoltage
- modifications to the appliance without the express approval of OPHARDT Hygiene
- use for any purpose other than that mentioned in the application

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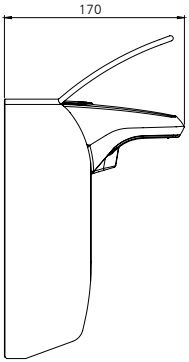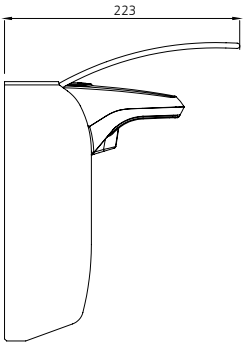
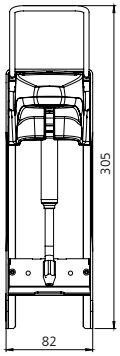
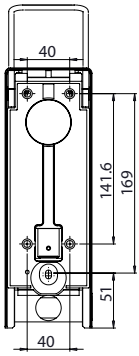
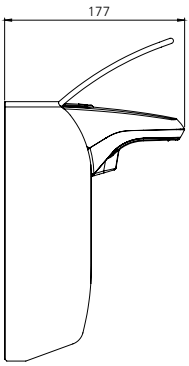
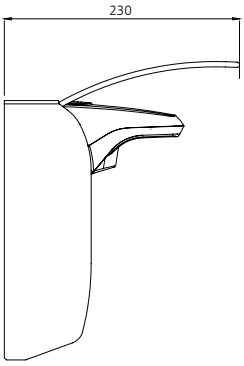
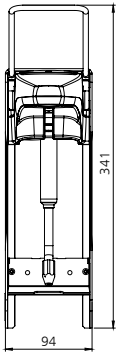
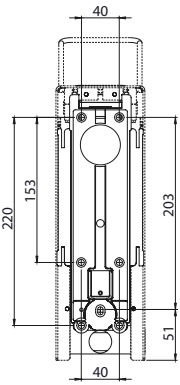
Ⓔ Important!

Please carefully read the operating instructions and save them for reference.
The proper function of this dispenser system can only be guaranteed if the procedures in these instructions are followed.

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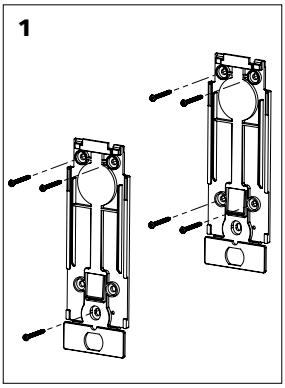
Ⓔ 1. Preparation

When using as a soap dispenser, install the dispenser as close as possible to the sink to prevent wet hands from dripping onto the floor. Wall plates already in use from the ingo-man® classic/plus product family may continue to be used. If mounting is necessary, please test wall for suitability before mounting.

Required space above the dispenser: 50-70 mm for the lever, and for dispenser removal and placement, as well as an additional 100 mm for operation. The required room below the dispenser, necessary for bottle removal, is about 130 for the IMS E/ELS and 170 mm for the IMS T/TLS. Optimal mounting height for use: 130-135 cm, measured from the floor to the upper edge of the dispenser housing.

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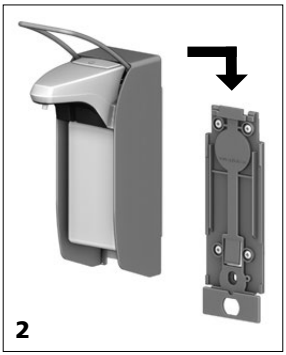
Ⓔ 2. Mounting of the dispenser

First remove the bottle: to do this tilt the bottom of the bottle forwards and then pull the bottle out of the dispenser. Remove the dispenser from the wall plate: press the latch on the back wall of the dispenser and slide the dispenser upwards, off the plate. Using at least 3 long screws provided, or other suitable fixing screws, attach the wall plate to the wall (see image 1). Slide the dispenser onto the wall plate from the top until the latch clicks into place (see image 2). If desired, the dispenser can be further secured by screwing the short screw provided through the hole in the dispenser's base bottle holder (shelf) into the backplate. Note: The dispenser cannot then be removed from the backplate without first removing this screw.

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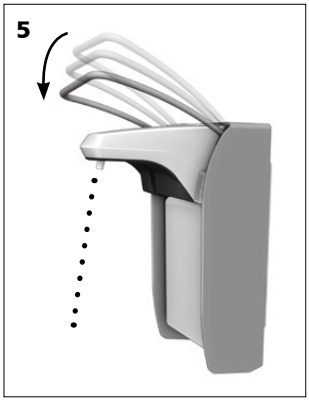(E) 3. Flexible Sealing Cap

The included flexible sealing cap ensures optimal sealing for various bottle sizes and prevents the loss of alcohol from the disinfectant solution. When fitting a pump to the dispenser, always ensure that the bottle sealing cap is pulled down the pump cylinder at least 1cm.

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(E) 4. Bottle Placement

Insert the filled bottle into the dispenser from the bottom, until snaps into place. Please also note point 15 on page 20.

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(E) 5. Initial Operation


For first use hold your hand under the nozzle and operate the lever several times until solution is dispensed.

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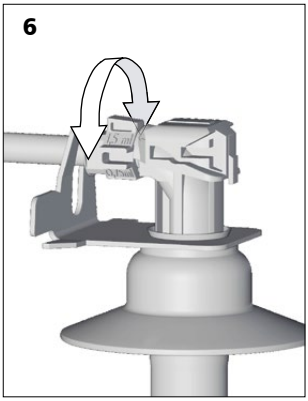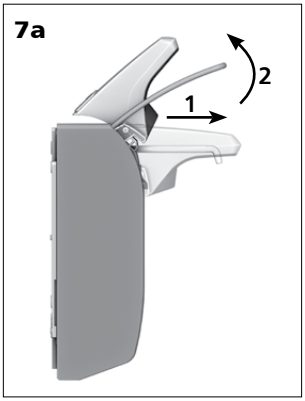
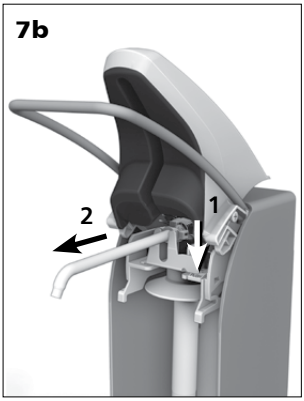
E 6. Dosage

The dispenser is pre-set to dispense the maximum dosage of either 1.5 ml or 1.0 ml/stroke, depending on the pump installed.

Important! Do not take apart the housing! Use only original replacement parts issued from the manufacturer. Only for the dosage of suitable solutions such as soaps, lotions and disinfectants. Caution notice! The use of unsuitable solutions can influence the function of the dispenser and may cause the solution to spray into eyes.

To change the dosage, first remove the pump from the dispenser (see step 7).
The dosage can be changed by turning the adjustable wheel. On first use, the wheel may be difficult to turn; if necessary, use a pair of flat-nose pliers. The set dosage can be read from the wheel, when looking at it from the front, right side. The upper indication of the dosage volume is relevant (see image 6).

Dosiermöglichkeiten je Pumpe - Dosage options depending on pump

Pumpe - Pump	(I)	(II)	(III)	(IIII)
1.5 ml (24, 25, 28)	1.5 ml	1.2 ml	1.0 ml	0.75 ml
K	1.5 ml	1.0 ml	0.7 ml	x
KSP	1.0 ml	x	x	x

Pumpe - Pump	(I)	(II)	(III)	(IIII)
DHP	1.5 ml	x	x	x
F	0.75 ml	x	x	x

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E 7. Pump Removal / Exchange

Gently pull face shield forward and then lift upwards. Next, press down on the pump lock with a finger, and pull the pump forwards (see images 7a + 7b).

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8

E 8. Cleaning / Sterilization

Remove the dispenser from the wall plate: to do this, press the latch on the back wall of the dispenser and slide the dispenser upwards, off the plate.

The ingo-man® smart dispenser (without SmartNose) and stainless-steel pumps (outlet pipe made from stainless steel), are autoclavable with superheated steam (134 °C, 2 bar, max. 20 min dwell time).

⚠ **Important!** SmartNose, wall plate, bottle and plastic pumps (K pump) are not autoclavable! For pump removal, see step 7. After autoclaving dispenser and the stainless steel pump, allow it to cool for at least an hour prior to next use.

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
E 9. Information for General Care

The Robert Koch-Institute recommends: “For both dispenser types (note: hand disinfectant and soap dispenser), the fixed outer and inner parts must be able to be reprocessed. This is also possible by machine. [...] Visible contamination of the outer housing and all accessible parts without further manipulation as well as soilings at the outlet must be removed by wipe disinfection. Daily wipe disinfection is recommended for operating levers. The scope and frequency of checks and the reprocessing of dispensers should be specified in the internal hygiene plan, e.g. every six months.” [Source: Robert Koch Institute, recommendations for hand hygiene, 2016]

In order to guarantee proper functioning of the dispenser regular manual cleaning or processing in cleaning-disinfection devices is necessary.

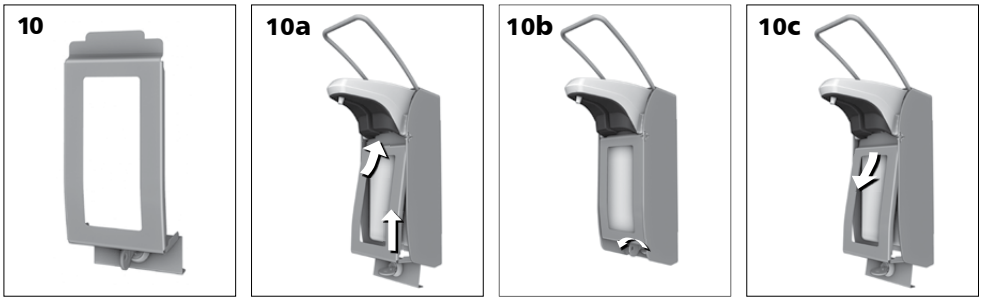Manual cleaning: Regularly remove adherent product residual. Clean the dispenser housing and the pump with hot water. Rinsing the pump with hot water. Disinfection of shroud is possible with appropriate surface disinfectants.
Processing in cleaning-disinfection devices: choose program „Dispenser cleaning" and use pH-neutral cleaning solutions. Important: Cleaning of plastic pumps up to max. 65 °C!
Processing in autoklave: the dispenser and stainless steel pump are autoclavable with superheated steam (134 °C, 2 bar, max. 20 min dwell time).
SmartNose, plastic wall plate, bottle, plastic pumps (K pumps) are not autoclavable.

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E 10. Locking plate (optional)

10a Insert the top of the locking plate into the slot under the face shield. The key should be in open position. Push the bottom of the locking plate up and forward to close it in place. Turn the key clockwise until the key is vertical to lock the dispenser. Retract the key from the lock

10b To open the locking plate insert the key and turn fully anticlockwise, by applying upward pressure on the locking plate, pull forward.

10c After the locking plate is disengaged from a backplate pull it down to detach it from the dispenser.

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
E 11. Scope of Delivery

Dispenser for soap and disinfectant
Fixing set: screws and plugs
Instruction manual

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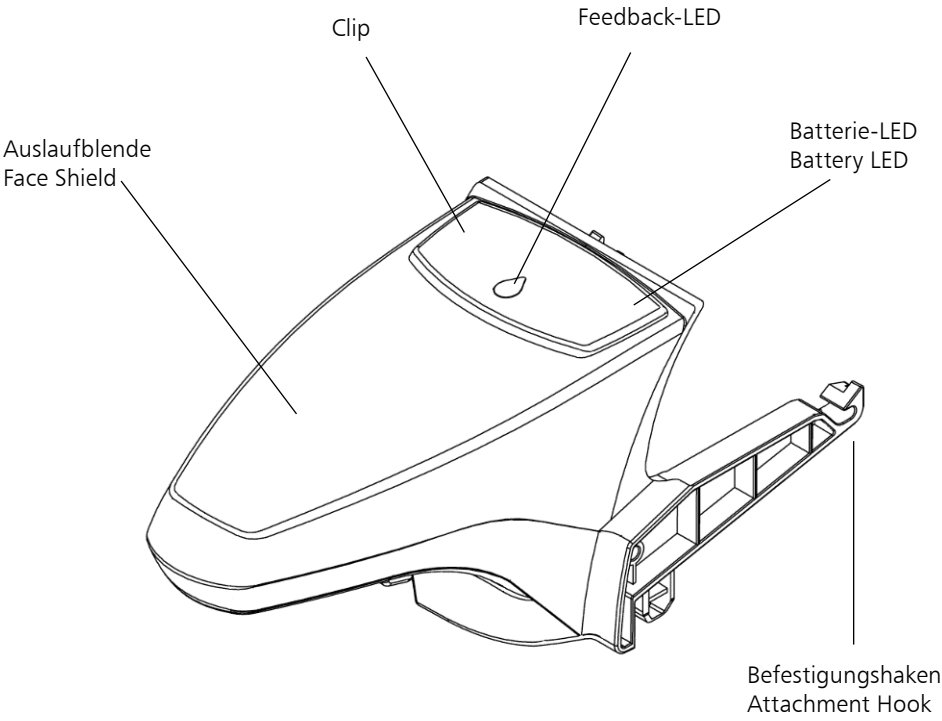E 12. Technical description:

The SmartNose automatically recognises the actuation of the dispenser and forwards information about the dosed quantity, the usage time and other sensor values via WLAN to a higher-level evaluation system.

When enough liquid has been dispensed, the feedback LED lights up green and signals to the user that enough liquid has been dispensed.

The unit uses a sensor to detect when a bottle is removed or inserted and signals this with a red or orange LED.

Device and WLAN settings can be configured using the NFC interface via the Kanary app.

Further processing of the data paired with a detailed overview, such as the calculation of the fill level or values for compliance, is possible with the help of the Kanary software. The software is not included in the scope of delivery.
For more information, please go to https://www.kanary.io or simply scan the QR code.

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13. Assembly

Installation preparation

During installation, make sure that no radiation-emitting object (e.g. X-ray machine, radio telephone, etc.) is in the immediate vicinity. Electromagnetic sources of interference and / or radiating electrical fields can impair the function.

Compatible with all IMP E/T dispensers.

The following requirements must be met in order to use all functions:

- A stable WiFi covering the area of use must be available in the IEEE 802.11 b/g/n standard
- The data from the SmartNose is collected in the Message Collector (OCG = Ophardt Communication Gateway). This can be installed locally or optionally provided by Ophardt hygiene in the cloud and transfers the data to the cloud platform Ophardt IoT.
- The SmartNose must be configured for use in the WLAN network.
- An internet-capable end device is required to evaluate the data (e.g. PC, laptop, tablet, smartphone).
- Please contact us if you have any questions about commissioning. We will be happy to help.

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Assembly

To install and commission the ingo-man SmartNose, please carry out the following steps in the order given:

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Open the outlet cover. To do this, pull slightly forwards and upwards.

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Lift the bottle slightly and remove it from the dispenser.

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
Hold down the pump lock with your thumb and pull out the pump to the front. Then swing the outlet cover down again.

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
To remove the outlet cover, pull the spring-loaded metal axle at the top of the dispenser forward with your index and middle fingers and unhook the fastening hooks by pivoting the rear of the outlet cover downward.

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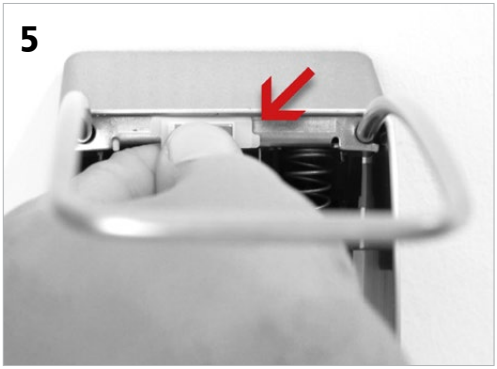Attach the magnet with the holder to the bracket base. To do this, press the magnetic holder into the opening provided and snap it into place.

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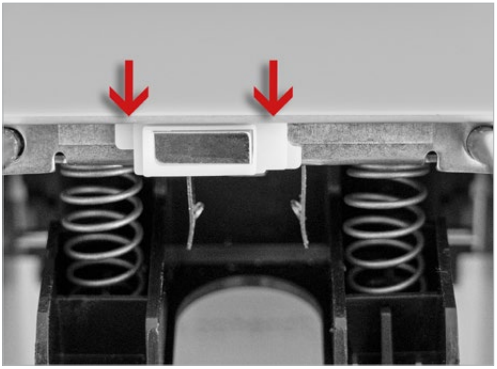Important: Ensure correct orientation of the magnet holder.

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Inserting the SmartNose: Pull the spring-loaded metal axle at the top of the dispenser forwards with your index and middle fingers and hook the fastening hooks of the outlet cover into place.

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Open the SmartNose. To do this, pull it slightly forwards and swing it upwards.

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
Insert the pump and close the SmartNose. Now check again that it is firmly seated. If the SmartNose is loose and can be moved, open it again and fasten it to the metal axle with both fastening hooks as described in point 7.

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Inserting the bottle: Insert the filled bottle into the dispenser from below and let it click into place. The dispenser is now ready for use.

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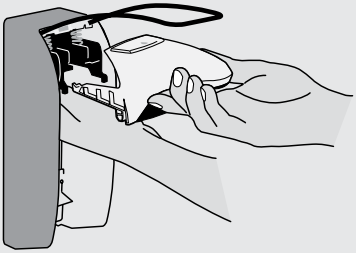Use appliance

Once the SmartNose has been installed on the dispenser, it adjusts to the environment and calibrates itself. During this process, the feedback LED flashes green. As soon as it goes out, the operating lever must be pressed all the way down 10 times. The SmartNose is then ready for operation.

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E 14. Make your disinfectant dispenser smart in just a few steps

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Install App

Scan QR-Code or search your app store.

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
Follow instructions

Follow the instructions shown in the app. After registration you have full access to the monitoring software KANARY.

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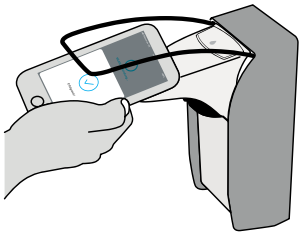
NFC Scan

Through the scan, the smart dispenser will be added to monitoring software KANARY.

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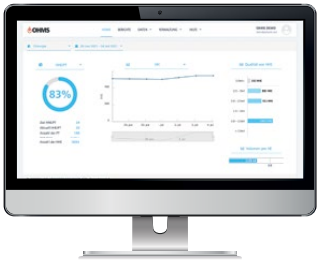E 15. Bottle exchange when using KANARY hygiene monitoring software

- If the locking plate is present it must first be removed.
- Lift up the empty bottle from below, draw it to the front and remove it from the dispenser.
- The Feedback-LED lights up red.
- Insert the filled-up bottle from below into the dispenser. Please pay attention to the correct position.
- The Feedback-LED lights up orange.

After a bottle change, the SmartNose attempts to establish a connection to the OCG (see page 14). If the connection is successful, all data from the SmartNose memory is transferred to the OCG. If the bottle is changed several times within 30 seconds, only one transmission takes place.

When changing bottles, it is assumed that a full bottle with the same content volume was used.

Ensure the correct configuration of the parameters “Bottle volume” and “Pump dosage” for an accurate calculation of the fill level.

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E 16. LED display

The feedback LED lights up green as soon as the required minimum amount for hand disinfection is reached within one hand hygiene event (standard: 3 ml). A hand hygiene event is the hand disinfection of one person with any number of actuations of the dispenser. A hand hygiene event is completed when a defined time has elapsed since the last actuation of the dispenser (default: 2 sec).

The feedback LED lights up green while the dispenser is calibrating after the SmartNose has been installed.

The feedback LED lights up orange when a bottle has been inserted into the dispenser.

The feedback LED lights up red when a bottle has been removed from the dispenser.

The red battery LED lights up during operation of the dispenser when the battery capacity is low.

When the SmartNose is removed from the dispenser, the green feedback LED and the red battery LED flash simultaneously. The flashing stops when the SmartNose is reinstalled or after approx. 1:30 min.

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E 17. Battery

If the battery LED lights up red during operation of the dispenser, the capacity is low.

Warning

Significant hazards may result from improper battery replacement. This unit contains a lithium battery that cannot be recharged.

- Never connect this battery to a charger.
- The battery must not be short-circuited.
- Dispose of empty batteries in line with local policy and protocol.

Dangerous goods

The lithium battery inside the SmartNose is a dangerous good. In case of return of dangerous goods, legal regulations apply. The return must be made in suitable packaging with full hazardous goods labelling. Do not ship units with damaged lithium batteries!

Spare battery

Operate the device only with an original battery, which is also available as a spare part.
(Lithium C Cell Battery 3.6 V, Art. No. 4402691 or Lithium battery ER 17505 3.6 V, Art. No. 3401863.
See page 27 for details.)

Battery change

Open the housing of the SmartNose. To do this, loosen the screws on both sides with a screwdriver.

Caution: The battery must not fall out of the housing.

Disconnect the battery from the circuit board. To do this, hold down the connector lock with your finger and carefully pull the connector off the circuit board.

Remove the empty battery from the housing.

Insert the new battery into the housing. The black cable on the battery must be routed between the two plastic tabs.

Snap the plug back into the connector on the PCB. The plug can only be inserted with the correct polarity.

Carefully close the SmartNose housing. The front tab of the lower part must be inserted between the two tabs of the upper part. Ensure that no area of the cable is pinched or squeezed.

Close the housing of the SmartNose. Use a screwdriver to tighten the screws on both sides.

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18. Troubleshoot

You will find some suggested solutions to possible faults in the adjacent table. Of course, you are also welcome to contact our service department: customercare@ophardt.com.

Properties	possible cause	Solution
The green feedback LED does not light up after a hand hygiene event.	The dispensed quantity is insufficient.	Activate the lever to dispense more fluid.
	Calibration is not completed/unit is not active.	Remove the SmartNose from the dispenser. Wait until the red battery LED and the green feedback LED flash. Reinstall the SmartNose.
	Battery is empty.	Insert new battery.
When inserting a container, the orange feedback LED does not light up.	Container was not recognised.	Remove the container and reinsert it.
	SmartNose (outlet cover) is open.	Close the SmartNose (outlet cover).
	Unit is defective.	Call customer service.
When removing the container, the red feedback LED does not light up.	SmartNose (outlet cover) is open.	Close the SmartNose (outlet cover).
	Unit is defective.	Call customer service.
The battery LED lights up red after each dispenser actuation.	The battery capacity is low.	Insert new battery.
	The battery has not been inserted correctly.	Insert the battery according to the instructions. Ensure correct polarity.
	The battery is unsuitable for the unit.	Only use original spare parts.
Usage data is not visible.	Please note that depending on the dispenser configuration, this must be pressed a few times before the transfer of the user data starts. It can take up to 30 minutes until the data of a completed transfer is visible in the application.	
	The wifi is not available.	Check the WLAN connection.
	Incorrect configuration of the WLAN or SmartNose.	Check configurations and correct if necessary.
The SmartNose cannot reach the message collector.	Message Collector is not running.	Start Message Collector
	The dispenser's server configuration is faulty.	Configure SmartNose.
	Message Collector has no connection to the Internet.	Establish an Internet connection.
	Message Collector proxy settings not adjusted.	Adjust proxy settings.
	A previous attempt to send has failed. The unit assumes that there is a faulty configuration or a network problem.	The SmartNose goes into a waiting phase (default: 8h). Only after this phase can a new transmission attempt start. To start an earlier transmission, remove the SmartNose from the dispenser, wait until the red battery LED and the green feedback LED flash. Then reinsert the SmartNose into the dispenser.

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19. Disposal

Packing

Put the packaging material in the appropriate recycling waste.

Device

When you finally take the unit out of service, please contact your recycling company or specialist dealer for information on disposal regulations.

Battery

The unit contains a lithium metal battery that must be disposed of properly.
Batteries must not be disposed of in household waste!

19. Technical data

Inside:	Board with sensors and WiFi module 802.11 b/g/n (WPA, WPA2 and WPA Enterprise), NFC module
Firmware:	Updates und Konfiguration „over-the-air“ (OTA)
Voltage:	3.6 V
Battery contained: in SmartNose E and SmartNose T	Lithium battery (Art. No. 4402691), 3.6 V, 9.0 Ah, transport code 3091, dangerous goods class 9, replaceable, not rechargeable, with cable and plug
Battery contained: in SmartNose INT E and SmartNose INT T	Lithium-Batterie (Art. No. 3401863), 3.6 V, 3.5 Ah, transport code 3091, dangerous goods class 9, replaceable, not rechargeable, with cable and plug
Protection class:	IP42
Conformity:	RoHS compliant
CE-marking:	EN 300328 Integration test WiFi ETSI EN 62311.2008 EMF ETSI EN 301489-1 EMV

NFC

OK FEEDBACK-LED

KANARY®

OPHARDT HYGIENE

OPHARDT HYGIENE-TECHNIK GmbH • Lindenau 27 • 47661 Issum • Germany
Phone +49 28 35 18-0 • Fax +49 28 35 18-77
Email customercare@ophardt.com • www.ophardt.com

4403088 · 26-01 · Rev. 2.0